


CH-100, Heating/Cooling Dry Block



Heat up and cool down times for CH-100

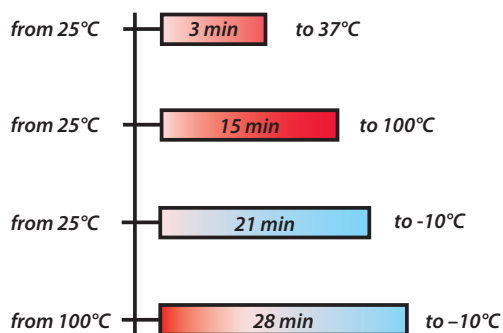

CH-100 is the result of combining two popular Biosan instruments:

1. Heating Dry block and
2. Cooling Dry block thermostat

The combined construction of aluminium block and Peltier element module cooled with the forced ventilation radiator provides fast changing of the cooling and heating modes.

CH-100 is a very effective instrument of sample preparation during enzyme reactions, hybridization reactions, DNA analysis.

Microprocessor controlled time and temperature. Simultaneous indication of set and actual temperature and time.

Specifications:

Temperature setting range	-10°C ... +100°C
Temperature control range	30°C below ambient ... +100°C
Setting resolution	0.1°C
Stability	±0.1°C
Digital timer with sound alarm	1 min – 96 hrs
Display	LCD, 2x16 signs
Overall dimensions	240 x 260 x 165 mm
Weight, not more	3.6 kg
External power supply	12V DC, 4.16 A

Blocks capacity:

Block CH-1	20x0.5 ml + 12x1.5 ml microtubes
Block CH-2	20x1.5 ml microtubes
Block CH-3	20x2 ml microtubes

Catalogue number:

CH-100 with block CH-1	BS-010410-BAI
CH-100 with block CH-2	BS-010410-CAI
CH-100 with block CH-3	BS-010410-UAI

Ice on block CH-2

